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INSTRUCTIONS FOR USE

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TECHNICAL NAME: Surgical Instrument

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TRADE NAME: Instrumental unarticulated, non-cutting, with the distal part in stainless steel, titanium and polyacetal (POM) - TRAUMEC.

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WARNING: Read carefully all of the instructions before use. Comply with all warnings and precautions mentioned in these instructions. The failure to comply with these recommendations may lead to complications.

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1. Product Identification

The surgical instruments are items specially developed for use in conjunction with implantable devices, assisting in installation or implantation. The devices of the family of articulated instrumental, not cutting, with distal part in stainless steel are invasive surgically transient use, reusable, designed and developed exclusively for use in implantable orthopedic devices used in surgical procedures.

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Carburo de tungsteno

2. Instrument Raw Materials

The following materials are employed on manufacturing the instruments:

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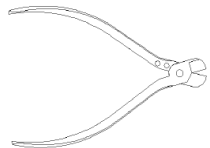
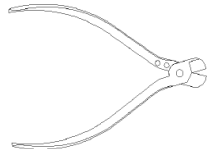
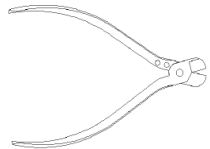
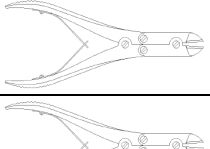
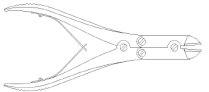
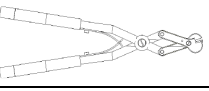
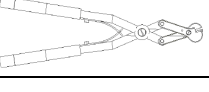
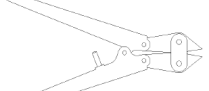
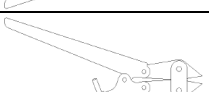
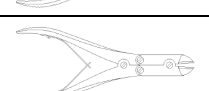
These Instructions also regard the following products:

CÓDIGO	DESCRIÇÃO	ILUSTRAÇÃO
PA.02.03.2756	Alicate de corte CPP Alicate de Corte CPP	
PA.02.03.2757	Alicate de corte COM Alicate de corte CPM	
PA.02.03.2758	Alicate de corte CPG Alicate de corte CPG	
PA.02.03.2759	Alicate de corte LPG Alicate de corte LPG	
PA.02.03.2760	Alicate de corte LPM Alicate de corte LPM	

PA.02.03.2761	Alicate de corte LPP Alicate de corte LPP	
PA.02.03.2762	Alicate de corte BPG Alicate de corte BPG	
PA.02.03.2763	Alicate de corte BPM Alicate de corte BPM	
PA.02.03.2764	Alicate de corte BPP Alicate de corte BPP	
PA.02.03.2765	Alicate de corte MPG Alicate de corte MPG	
PA.02.03.2766	Alicate de corte MPM Alicate de corte MPM	
PA.02.03.2767	Alicate de corte MPP Alicate de corte MPP	
PA.02.03.2768	Alicate de corte grande tipo tesoura Alicate de Corte Grande Tipo Tijera	
PA.02.03.2769	Alicate de corte médio tipo tesoura Alicate de Corte Medio Tipo Tijera	
PA.02.03.2770	Alicate de corte pequeno tipo tesoura Alicate de Corte Pequeño Tipo Tijera	
PA.02.03.2771	Alicate de corte dupla ação grande Alicate de Corte Doble Acción Grande	

PA.02.03.2772	Alicate de corte dupla ação médio <i>Alicate Ge Corte Doble Acción Mediano</i>	
PA.02.03.2773	Alicate de corte dupla ação pequeno <i>Alicate de Corte Doble Acción Pequeño</i>	
PA.02.03.2774	Alicate de corte rebarbador grande <i>Alicate de Corte Desbarbador Grande</i>	
PA.02.03.2775	Alicate de corte rebarbador médio <i>Alicate de Corte Desbarbador Mediano</i>	
PA.02.03.2776	Alicate de corte rebarbador pequeno <i>Alicate de Corte Desbarbador Pequeño</i>	
PA.02.03.2777	Alicate de corte tipo tesoura rebarbador pequeno <i>Alicate de Corte Tipo Tijera Desbarbador Pequeño</i>	
PA.02.03.2778	Alicate de corte tipo tesoura rebarbador médio <i>Alicate de Corte Tipo Tijera Desbarbador Mediano</i>	
PA.02.03.2779	Alicate de corte tipo tesoura rebarbador grande <i>Alicate de Corte Tipo Tijera Desbarbador Grande</i>	

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4. Operation principle

The Instruments has the working principle to assist the surgeon in surgical procedures.

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5. Directions

Preoperative: The instrument selection is an integral part of the surgical planning and must be conducted through a formal medical request indicating the intended technique, as well as characteristics of the implant to be used and specification of components which are part of the instrument system. It is of utmost importance to conduct a thorough inspection in each component, paying attention to the usage and cleaning conditions. Sterilization is mandatory and must present proven efficiency. The surgical instrumentation technique review before effective usage increases procedure efficiency.

Intraoperative: The instruments are exclusively intended to assist surgeons, and shall never be an integral part of implants, regarding the permanence inside the patient's body after the procedure. They must be selected and harmonized exclusively for the device to be implanted, with the purpose to ensure orthopedic implant adequacy on the implant location.

Instruments from different manufacturers may be used in association only when the compatibility and adequacy between them have been previously assessed.

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PROCEDURES IN ORDER TO USE AND REUSE THE MEDICAL PRODUCT

The application of surgical instruments must always be performed under technical guidance and restricted to clinical and hospital environments complying with the following safety measures:

Handling and transportation: The instruments must be transported and handled as to prevent any damages or alterations of their characteristics. They must be handled carefully, in small batches, as to prevent bumping or dropping. Any instruments that has been dropped or handled in an unsuitable manner, or even if it is suspected that it may have suffered any damages, must be isolated and submitted to the qualified technical supervisor of the institution in order to be inspected, even if it has already been through this step.

Inspections: Only the instruments that have previously undergone technical inspection may be used.

Technical Inspection: Before being made available for use, the instruments, including the set assembly, must undergo a technical inspection, performed by a qualified technical supervisor. The unapproved parts must be sorted for review and maintenance by supplier or disposed. The inspection must check whether the characteristics regarding preservation and operability of the instruments, including superficial aspects such as stains, oxidations and damages, besides characteristics regarding each instrument, such as ease of articulation, apprehension capacity, cutting capacity and edge alignment.

Sterilization: The instruments are supplied fully decontaminated, clean and non-sterile; they must be sterilized before use. TRAUMEC recommends the sterilization methods described below.

The parameters suitable for sterilization processes (physical or chemical) for each equipment and volume must be analyzed and conducted by people trained and specialized in sterilization processes, ensuring the complete efficiency of this procedure.

Cycle	Temperature	Exposure Time
Regular (pressure of 1 atm)	121°C (250°F)	30 minutes
Regular (1 atm pressure)	132°C (270°F)	15 minutes
Gravity	132°C (270°F)	45 minutes
High-Vacuum	132°C (270°F)	4 minutes

NOTE: Time must only be counted when the sterilization chamber reaches the desired temperature. The instruments must be sterilized in autoclaves before use, according to a validated method and in compliance with NBR ISO 17665 Part 1: 2010.

Reusing: The process for surgical instrument reusing involves, at least, five basic steps: pre-cleaning, decontamination, washing, rinsing and drying. It is recommended that each instrument is immediately cleaned after the surgical procedure it has been used in, preventing the hardening of residues resulting from this procedure. The cleaning must also be standardized, preventing the propagation of contaminants and instrument damage. Each manual cleaning procedure must be performed while using suitable personal protection equipment. When cleaning instruments in automated equipment, the manufacturer instructions must be strictly followed, especially regarding the quality of products and water involved in the cleaning procedure. The instruments, whenever relevant, must be inserted in the equipment while opened or disassembled.

Under no circumstances whatsoever wire brushes, steel wool or any other abrasive products, even scouring products, must be used in order to remove residues resulting from any steps of the cleaning process. Do not use aggressive cleaning agents such as mineral and acid agents (sulphuric, nitric).

It must be guaranteed that the instrument, as well as its components, whenever relevant, are free of any preserving agents, as well as from any residues that may result from storage or repairing procedure. The presence of products not soluble in water may result in the creation of physical barriers, protecting microorganisms from germicides, as well as it may allow unwanted residue retention after use. The water quality is a crucial factor both to the cleaning process and to the instrument preservation. The presence of particulates, the concentration of elements or chemical substances, and the pH unbalance may cause instrument deterioration during the cleaning process.

The combination of some of these parameters may lead to mineral precipitation deposits, which cannot be eliminated in the organic matter incrustation removal phase, as well as it may start the corrosion process of the stainless steel, as in the case of excessive chloride concentration. It is recommended that the water used to wash the instruments complies with the quality requirements established in the sterilization process.

Note: All instruments must be cleaned right after finishing the surgical procedure; by doing so, the hardening of liquids resultant from the surgical procedure is avoided. The entire cleaning process must be performed carefully, preventing drops and bumps that may compromise the instrument operability.

Pre-cleaning: The instrument must be preferably submerged and, whenever possible, opened or disassembled in a suitable container with water and a preferably enzymatic detergent at room temperature. Afterwards, it must be strictly washed under running and preferably lukewarm water. This phase must always be performed with water under 45°C, since elevated temperatures cause protein coagulation, turning the process of removing incrustations from the instrument a bit harder.

Decontamination: It is performed through submerging, whenever possible, an open or disassembled instrument, in a suitable container, with a solution of disinfectant in water at room temperature (chemical disinfection), or in a warmed bath (thermochemical disinfection). The amount of time in which the instrument must be submerged depends on the operation temperature, dilution and the type of disinfectant used.

Washing: The pieces must be thoroughly brushed, with a soft bristle brush, paying special attention to articulations, grooves and gear racks. The instrument, whenever relevant, must be disassembled and each component must be washed separately. Pay special attention to areas that are difficult to reach, where organic tissue, secretions or disinfection solutions may be retained or deposited.

Rinsing: The instrument must be plentifully rinsed under running water; articulated instruments must be opened and closed several times during the rinsing process. It is recommended to use warm water in order to rinse the instrument.

Drying: You must ensure that the drying processes do not introduce humidity, particles or fuzz over instrument surface. Pay special attention to articulations, grooves and gear racks. It is recommended that the tissue is absorbent, soft and that each part of a disassembling instrument is dried separately; and if the instrument presents cavities or slots, their inside must be completely dried.

Disposal: The disposal of disqualified parts must be performed under technical assessment and guidance. After the replacement, destroy all damaged components, preventing further improper use.

When there is the need to dispose of the instrument, it must be immediately disabled, in order to prevent its inadvertent use.

The disposal of instruments must comply with the standards regarding disposal of contaminant medical waste.

We recommend that the parts are filed, bent or cut in order to be disabled. For the disposal of instruments, follow your local legislation in force regarding the discard of potentially contaminant products.

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6. Storage conditions

The instruments must be stored in a clean, dry and ventilated place, at room temperature and away from direct light. The special storage, handling and preservation conditions of the product must be observed, in order to guarantee that the components remain intact until used in the surgical procedure. Safety measures regarding receipt, storage, transport, cleaning and batch references preservation must be adopted along the best practices of storage and distribution of medical products.

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7. Handling and transportation conditions

The instruments must be transported and handled as to prevent any damages or alterations of their characteristics. The instruments must be managed and handled carefully, in small batches, preventing bumping or dropping. Any instruments that has been dropped or handled in an unsuitable manner or even if it is suspected that it may have suffered any damages must be isolated and submitted to the responsible technician, qualified for performing the inspection, even if it has already been through this step.

Note: Any products that has been dropped or unsuitably managed, or even if it is suspected it may have been suffered any damages ,must be sorted and segregated.

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9. Warnings

There are surgical instruments suitable to every step of a surgery. The usual wear, exertion of excessive forces and use of these instruments for ends which do not befit the usage indications may compromise the progress of surgical procedure and damage the implant. Also, using instrument supplied by different manufacturers may present risks, such as unsuitable fixation of the implant and instrument damage. The instruments are metallic components subject to important mechanical requirements during repeated use

for variable and indefinite periods of time; therefore, the inspection and review of the usage conditions of instruments or their parts is required. In the event there is a varying performance, precision loss, instability or lack of sharpness, the part must be immediately replaced.

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10. Precautions

The instruments must be kept in their original packages until sterilized and used. Clean it correctly after every use, in order to prevent incrustations or corrosions.

Only professionals who are specialized and trained on the correspondent surgical technique are able to use these instruments.

At every use, check whether the instruments have suffered any damages.

Always use the instruments suitable for every type of implant; do not make attempts to substitute any elements for others which are not appropriate for the intended use.

The instruments must not be stored along with chemical products that may produce corrosive vapors, causing eventual damage to them.

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www.traumec.com.br

Rev.: 00

Technical name: SURGICAL INSTRUMENTS

Trade Name: INSTRUMENTAL UNARTICULATED, NON-CUTTING, WITH THE DISTAL PART IN STAINLESS STEEL, TITANIUM AND POLYACETAL (POM) - TRAUMEC.

Business Model: PLIER CUTTER CPP

Code: PA.02.03.2746

Lot: XXXXXX

Quantity: 01 pc

ANVISA Registration No: 80455630045

Manufacture: XX / XX / XXXX

Validity: XX / XX / XXXX

Material: XXXXXX

Technical Director: José Luiz Caritá - CREA-SP - 0685038754

PRODUCT NOT STERILE

SPECIAL CONDITIONS OF STORAGE, STORAGE, HANDLING MEDICAL PRODUCT: SEE INSTRUCTIONS.

ANVISA Registration No 80455630045

Technical Director: José Luiz Caritá - CREA-SP - 0685038754

Manufactured by: Traumec Orthopedic Implants and Technology Import & Export Ltda.

Address: 1A JC Street, 138, Subdivision: Centennial Garden

City: Rio Claro | UF: SP | CEP -13503-510

CNPJ: 09.123.223/0001-10 - Brazilian Industry

Customer Service - Phone / Fax: +55 19 3522-1177

E-mail: sac@traumec.com.br

www.traumec.com.br

Rev.: 00